


# Health N Massage®

## NECK MASSAGER WITH NATURE'S SOUNDS

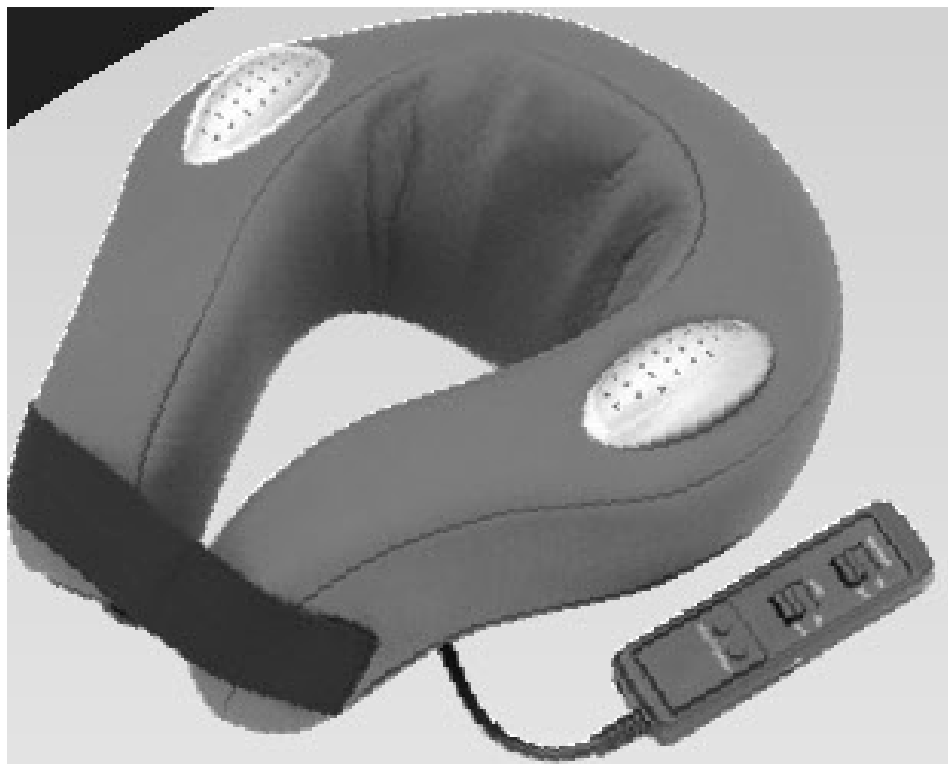

**ITEM NO: 900806**

**MODEL NO: KM 7600AA**

## **IMPORTANT SAFETY INSTRUCTIONS**

When using an electrical appliance, basic precautions should always be followed, including the following:

Read all instructions before using this appliance.

Keep this appliance out of the reach of children and infirm persons. Young children should be supervised to ensure that they do not play with the appliance.

**DANGER** – To reduce the risk of electric shock:

1. Always unplug this appliance from the electrical outlet immediately after using and before cleaning.
2. DO NOT reach for a product that has fallen into water. Unplug immediately.
3. DO NOT use while bathing or in the shower.
4. DO NOT place or store appliance where it can fall or be pulled into a bath or sink. DO NOT place in or drop into water or other liquid.

**WARNING** – To reduce the risk of burns, fire, electric shock, or injury to persons:

1. An appliance should never be left unattended when plugged in. Unplug from mains when not in use and before putting on or taking off parts.
2. DO NOT operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury.
3. Use this appliance only for its intended use as described in this manual, DO NOT use attachments not recommended by the manufacturer.
4. DO NOT carry this appliance by supply cord or use cord as a handle.
5. Keep the cord away from heated surfaces.
6. DO NOT use outdoors.
7. Never drop or insert any object into any opening.
8. To disconnect, turn all controls to the off position, then remove plug from mains.


## **PRODUCT CARE AND MAINTENANCE**

- Place massager in a safe, damp cloth. DO NOT use abrasive cleaner.
- Never submerge the unit in any liquid. Keep away from all solvents and harsh detergents.
- DO NOT dry clean.
- DO NOT attempt to repair this massager. There are no user-serviceable parts.

## **HOW TO USE**

At home or at the office:

1. Place unit comfortably around neck.  
Adjust the neck massager to the part of the neck you wish to treat. Attach the self holding strap across the front of the unit. (see Figure A)
2. Turn on the massage motor by sliding the massage control button to the "ON" (I) position. (See Figure B)
3. Relaxing heat can be applied by sliding the heat control button to the "ON" (I) position. The heat setting operates only when the unit is powered the AC adaptor.
4. The unit provides magnetic therapy with magnets which are located within the massager. Magnetic therapy helps improve blood flow and promotes the body's natural healing process.
5. Eight soothing, natural sounds are generated from two build-in speakers. Sounds can be as loud or as soft as required by using the volume control dial located on the side of the handset.
6. To activate sounds, push the ON/OFF (I/O) button on the control panel, then push the SELECT & button to choose from sounds including song birds, rural sounds, white noise, ocean wave, summer night, brook, running stream, rainforest. (See Figure D)
7. DO NOT use for more than 15 minutes continuous use.
8. Unit can either run by AC adaptor or batteries. Four "AA" size batteries are required for battery operation.

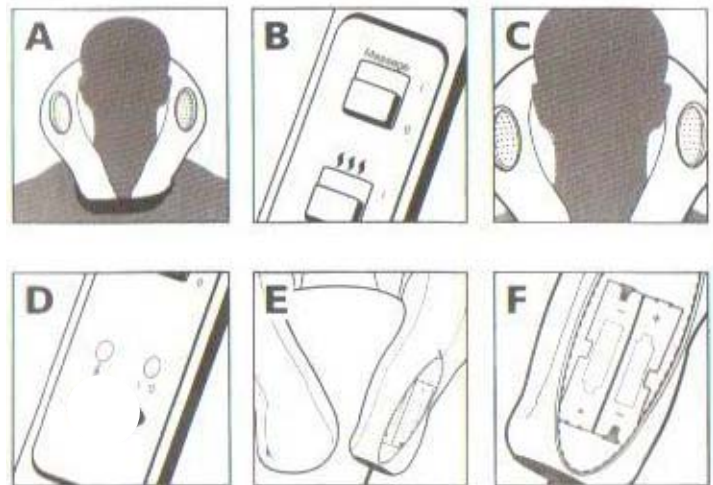

## **BATTERY INSTALLATION**

Four "AA" size batteries are required for operation

1. To install batteries, slide open one side zip on the back of the unit. (See Figure E)
2. Insert two batteries in the plastic battery compartment as indicated. (see Figure F)
3. Close the zip.
4. Repeat this procedure for the other side of the unit.

## **GUARANTEE**

### **WHAT'S COVERED**

**H. E. Group Ltd** guarantees **HEALTH N MASSAGE – Neck Massager with Nature's Sounds (Item no: 900806)** for a full 12 months replacement period from the date of purchase, which is only valid with the presentation of the original proof of purchase receipt. This guarantee covers materials and workmanship defects only.

### **WHAT IS NOT COVERED**

This guarantee does not apply to defects arising from abuse, alteration or unreasonable use of the product. In addition, it does not cover general deterioration due to usage.

This guarantee does not apply to the **accessories**, due to general deterioration.

### **WHAT H. E. GROUP LTD WILL DO**

**H. E. Group Ltd** will replace the product, if a defect is discovered whilst under the guaranteed period, subject to verification of defect or malfunction, date and proof of purchase.

### **CUSTOMER CARE**

*All enquiries may be directed to our toll free customer service number:  
1800 079 124 or via email [info@hegroup.com.au](mailto:info@hegroup.com.au)*